

APPENDIX B

Molecular Probes™

Invitrogen detection technologies

Product Information

Revised: 22-August-2005

Zenon® Mouse IgG Labeling Kits

Quick Facts

Storage upon receipt:

- 2–6°C
- Protect from light

Ex/Em: See Table 1

Introduction

The Zenon® Mouse IgG Labeling Kits provide a fast, versatile and reliable method for producing antibody conjugates, even with very small (submicrogram) amounts of starting material. Antibody conjugates formed using Zenon® technology may be used to stain cells in any application where a directly labeled primary antibody is suitable, including flow cytometry, imaging, high throughput and other applications. Moreover, this technology simplifies applications that previously were time consuming or not practical, such as the use of multiple mouse-derived antibodies in the same staining protocol.

Each of the Zenon® Mouse IgG Labeling Kits is designed for use with a particular mouse monoclonal antibody isotype: IgG₁, IgG_{2a} or IgG_{2b}. The available labels consist of a large selection of premium fluorescent dyes, including our entire range of Alexa Fluor® dyes, as well as conventional dyes, R-phycoerythrin (R-PE), allophycocyanin (APC), Alexa Fluor® dye-R-PE or Alexa Fluor® dye-APC tandem constructs, enzymes and biotin (Table 1). We also offer Zenon® Tricolor Mouse IgG Labeling Kits, each with a selection of three different Zenon® mouse IgG labeling reagents (Table 2).

Zenon® labeling technology utilizes a fluorophore-, biotin- or enzyme-labeled Fab fragment directed against the Fc portion of an intact IgG primary antibody in order to form a labeling complex (Figure 1). The labeled Fab fragments have been affinity purified during their preparation to ensure their high affinity and selectivity for the Fc portion of the primary antibody. Because this labeling is based on immunoselectivity, the Zenon® labeling method does not require the removal of exogenous proteins such as serum albumin or amine-containing buffers from the antibody prior to forming the complex. Crossreactivity is low with antibodies from other species.

Formation of the Fab-antibody complex occurs in less than 5 minutes, and nearly all of the primary antibodies in the mixture are labeled (Figure 2). Complexes formed using this technology

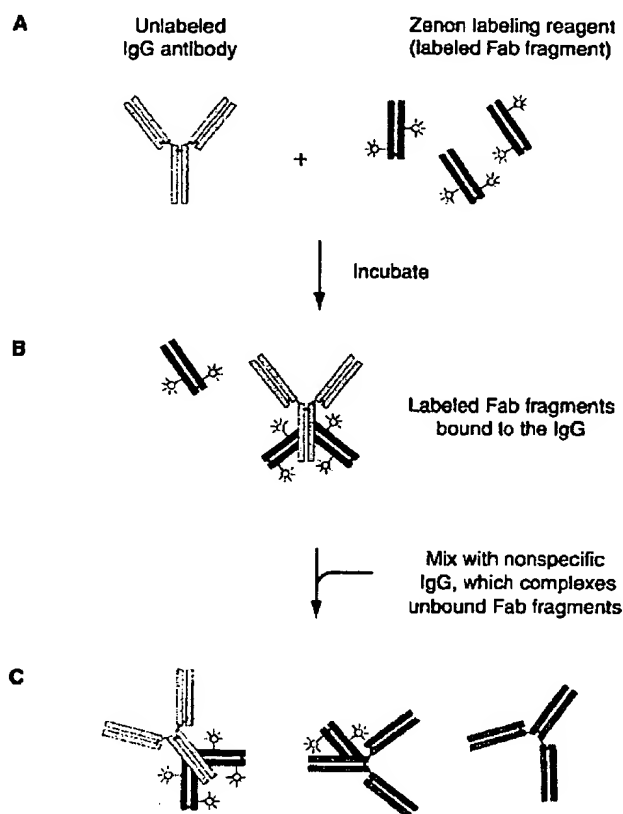


Figure 1. The Zenon® labeling scheme. An unlabeled IgG is incubated with the Zenon® labeling reagent, which contains a fluorophore-labeled Fab fragment (A). The labeled Fab fragment binds to the Fc portion of the IgG antibody (B), and excess Fab fragment is neutralized by the addition of a nonspecific IgG (C). The addition of nonspecific IgG prevents cross-labeling of the Fab fragment in experiments where multiple primary antibodies of the same type are present. Note that the Fab fragment used for labeling need not be coupled to a fluorophore, but could instead be coupled to an enzyme or to biotin.

display fluorescence intensity or enzymatic activity similar to that of directly labeled primary antibodies. In addition, the extent of antibody labeling (and thus the fluorescence intensity or enzymatic activity of the probe) can be adjusted by varying the amount of Zenon® labeling reagent that is added, i.e. by varying the molar ratio of labeled Fab fragment to primary antibody.

Molecular Probes maintains an e-mail address (probeszenon@invitrogen.com) specifically dedicated to the Zenon® labeling technology. For technical questions and/or comments about these products, please send an e-mail to the above address, or contact our Technical Service department.

MP 25000

Zenon® Mouse IgG Labeling Kits

Materials

The Zenon® Mouse IgG Labeling Kits contain sufficient reagents for either 50, 25 or 10 labelings, depending upon the kit (see *Contents*, below). One labeling is defined as the amount of Zenon® labeling reagent required to label 1 µg of an intact, affinity-purified mouse IgG antibody at a Fab:antibody molar ratio of 3:1 (note A).

Contents

Zenon® Mouse IgG Labeling Kits for labeling with a low molecular weight dye or biotin:

- Zenon® mouse IgG labeling reagent (Component A), 250 µL
- Zenon® blocking reagent (mouse IgG) (Component B), 250 µL

Sufficient reagents are provided for 50 labelings.

Zenon® Mouse IgG Labeling Kits for labeling with R-phycoerythrin (R-PE), allophycocyanin (APC) or an enzyme:

- Zenon® mouse IgG labeling reagent (Component A), 125 µL
- Zenon® blocking reagent (mouse IgG) (Component B), 125 µL

Sufficient reagents are provided for 25 labelings.

Zenon® Mouse IgG Labeling Kits for labeling with Alexa Fluor® dye-R-PE or Alexa Fluor® dye-APC tandem constructs:

- Zenon® mouse IgG labeling reagent (Component A), 50 µL
- Zenon® blocking reagent (mouse IgG) (Component B), 50 µL

Sufficient reagents are provided for 10 labelings.

Zenon® Tricolor Mouse IgG Labeling Kits (see Table 2):

- Zenon® mouse IgG labeling reagent of one fluorescent color (Component A), 50 µL
- Zenon® mouse IgG labeling reagent of a second fluorescent color (Component B), 50 µL
- Zenon® mouse IgG labeling reagent of a third fluorescent color (Component C), 50 µL
- Zenon® blocking reagent (mouse IgG) (Component D), 150 µL

Sufficient reagents are provided for 10 labelings with each of the three included labeling reagents.

Zenon® mouse IgG labeling reagents are labeled goat Fab fragments selective for the Fc portion of mouse IgG antibodies of the indicated IgG isotype—IgG₁, IgG_{2a} or IgG_{2b}—depending on the kit. Zenon® mouse IgG labeling reagents conjugated to a low molecular weight dye or biotin are supplied at 200 µg Fab fragment/mL in 0.1 M sodium phosphate, 0.1 M NaCl, pH 7.5, containing 5 mM sodium azide. Zenon® mouse IgG labeling reagents conjugated to R-phycoerythrin, allophycocyanin, Alexa Fluor® dye-R-PE or Alexa Fluor® dye-APC tandem constructs or alkaline phosphatase are supplied at 200 µg Fab fragment/mL in 0.1 M sodium phosphate, 0.1 M NaCl, pH 6.8, containing 5 mM sodium azide. Zenon® mouse IgG labeling reagents

Table 1. Molecular Probes' Zenon® Mouse IgG Labeling Kits.

Label	Ex/Em *	Catalog Number		
		Mouse IgG ₁	Mouse IgG _{2a}	Mouse IgG _{2b}
Alexa Fluor® dyes				
Alexa Fluor® 350	346/442	Z25000	Z25100	Z25200
Alexa Fluor® 405	402/421	Z25013	Z25113	Z25213
Alexa Fluor® 430	434/539	Z25001		
Alexa Fluor® 488	495/519	Z25002	Z25102	Z25202
Alexa Fluor® 532	531/554	Z25003		
Alexa Fluor® 546	558/573	Z25004	Z25104	Z25204
Alexa Fluor® 555	555/565	Z25005	Z25105	Z25205
Alexa Fluor® 568	578/603	Z25006	Z25106	Z25206
Alexa Fluor® 594	590/617	Z25007	Z25107	Z25207
Alexa Fluor® 647	650/668	Z25008	Z25108	Z25208
Alexa Fluor® 660	663/690	Z25009	Z25109	Z25209
Alexa Fluor® 680	679/702	Z25010	Z25110	Z25210
Alexa Fluor® 700	696/719	Z25011		
Alexa Fluor® 750	752/779			
Conventional dyes				
Marina Blue®	365/460	Z25040		
Pacific Blue™	410/455	Z25041	Z25156	
Fluorescein	494/518	Z25042		
Oregon Green® 488	496/524	Z25043		
Texas Red®-X	595/615	Z25045		
Biotins				
Biotin-XX	NA	Z25052	Z25152	Z25252
DSB-X™ † biotin	NA	Z25053		
R-phycoerythrin, allophycocyanin and tandem constructs				
R-phycoerythrin (R-PE)	496 ‡/578	Z25055	Z25155	Z25255
Alexa Fluor® 610-R-PE	496 ‡/630	Z25020		
Alexa Fluor® 647-R-PE	496 ‡/668	Z25021	Z25121	Z25221
Alexa Fluor® 680-R-PE	496 ‡/702	Z25022		
Allophycocyanin (APC)	650/660	Z25051	Z25151	Z25251
Alexa Fluor® 700-APC	650/723	Z25030		
Alexa Fluor® 750-APC	650/775	Z25031		
Enzymes				
Horseradish peroxidase	NA	Z25054	Z25154	Z25254
Alkaline phosphatase	NA	Z25050	Z25150	Z25250

* Approximate fluorescence excitation and emission maxima, in nm. † DSB-X is our trademark for desthiobiotin-X. ‡ Additional absorption peaks are present at 546 and 565 nm. NA = not applicable.

* Approximate fluorescence excitation and emission maxima, in nm. † DSB-X is our trademark for desthiobiotin-X. ‡ Additional absorption peaks are present at 546 and 565 nm. NA = not applicable.

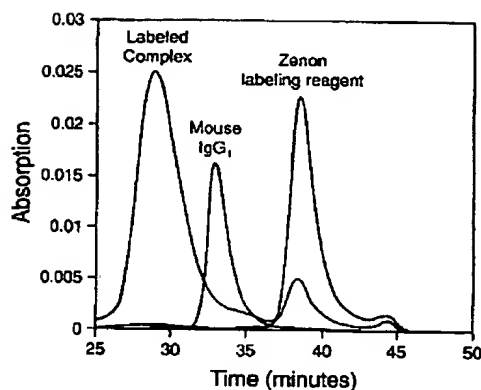


Figure 2. High-performance size-exclusion chromatographic analysis of the Alexa Fluor® 488 Zenon® labeling reagent binding to a mouse IgG₁ antibody. The Zenon® labeling reagent, alone, appears as a peak at 38 minutes; the mouse IgG₁, alone, at 33 minutes. When mixed together at a molar ratio of ~5:1 (Zenon® labeling reagent to IgG₁), the IgG₁ is quantitatively converted to a labeled complex that appears as a peak at 29 minutes.

conjugated to horseradish peroxidase are supplied at 200 µg Fab fragment/mL in 0.1 M sodium phosphate, 0.1 M NaCl, pH 6.8, containing 0.02% thimerosal. The Zenon® blocking reagent (mouse IgG) is supplied at 5 mg/mL in phosphate-buffered saline, pH 7.2, containing 5 mM sodium azide.

Storage and Handling

Upon receipt, store all Zenon® Mouse IgG Labeling Kits at 2–6°C, and protect fluorescent components from light. For long-term storage, Zenon® mouse IgG labeling reagents that contain a low molecular weight fluorophores or biotin and the Zenon® blocking reagent can be divided into single-use aliquots and frozen at ≤–20°C. DO NOT FREEZE Zenon® mouse IgG labeling reagents containing R-phycoerythrin, allophycocyanin, Alexa Fluor® dye-R-PE or Alexa Fluor® dye-APC tandem constructs, horseradish peroxidase or alkaline phosphatase.

Zenon® Complex Formation

The following protocol is for labeling 1 µg of antibody with a Zenon® mouse IgG labeling reagent to obtain a 3:1 molar ratio of Fab to antibody target (note A). This molar ratio is a suggested starting point and represents the minimum ratio for adequate labeling in most applications; individual experiments may require higher molar ratios in order to obtain satisfactory signal (note B). For larger or smaller quantities of antibody, the amounts of the reagents specified in this protocol can be scaled accordingly. The Zenon® mouse IgG labeling reaction does not require the removal of bovine serum albumin (BSA) or other stabilizing proteins that may be present in antibody preparations. Antibodies contained within ascites fluid or hybridoma supernatants can also be directly labeled and do not require purification of the antibody either prior to or after labeling.

NOTE: Before beginning the labeling protocol, verify that the isotype of the mouse antibody matches that of the Zenon® Mouse IgG Labeling Kit. The Zenon® mouse IgG labeling reagents in the kits are isotype specific and are not recommended for labeling antibodies that are not of the corresponding mouse IgG isotype.

1.1 Prepare 1 µg of antibody in a suitable buffer, such as phosphate-buffered saline (PBS). The volume is not crucial, provided it is ≤20 µL. If using an antibody where the concentration is not specified, refer to note C. If using an impure antibody preparation, e.g. ascites fluid or hybridoma supernatant, refer to note D.

1.2 Add 5 µL of the Zenon® mouse IgG labeling reagent (Component A) to the antibody solution. (If using a Zenon® Tricolor Kit, use the desired labeling reagent from Component A, B, or C.)

1.3 Incubate the mixture for 5 minutes at room temperature.

1.4 Add 5 µL of the Zenon® blocking reagent (Component B, or Component D if using a Zenon® Tricolor Kit) to the reaction mixture. The blocking reagent is not isotype specific and is the same in all Zenon® Mouse IgG Labeling Kits. If the application uses only a single labeled antibody, refer to note E.

1.5 Incubate the solution for 5 minutes at room temperature.

The complexes are now ready and should be applied to samples within approximately 30 minutes. If longer storage is required, refer to note F.

General Application Tips

Conjugate utility. Mouse IgG antibodies labeled using Zenon® technology are expected to be suitable in all applications where a directly labeled antibody can be used. Multiple mouse IgG antibodies labeled with Zenon® reagents can be used in one experiment either sequentially or as a single staining mixture.

Conjugate stability. Once the conjugates have been formed and excess Fab taken up by the blocking reagent, the labeled complexes should be used within approximately 30 minutes.

Working concentration. For applications where a directly labeled primary antibody is typically used, the antibody labeled using Zenon® technology can generally be used at a similar or higher concentration (1.5–3-fold).

Table 2. Molecular Probes' Zenon® Tricolor Mouse IgG Labeling Kits.

Tricolor Labeling Kit	Labels *	Catalog Number†		
		Mouse IgG ₁	Mouse IgG _{2a}	Mouse IgG _{2b}
Kit #1 (for imaging)	Alexa Fluor® 488 Alexa Fluor® 555 Alexa Fluor® 647	Z25060	Z25160	Z25260
Kit #2 (for imaging)	Alexa Fluor® 350 Alexa Fluor® 488 Alexa Fluor® 594	Z25070	Z25170	Z25270
Kit #3 (for flow cytometry)	Alexa Fluor® 488 R-phycoerythrin (R-PE) Alexa Fluor® 647-R-PE	Z25080	Z25180	Z25280

* See Table 1 for fluorescence excitation and emission maxima.

Low or no signal. If suitable controls verify that the primary antibody is binding to the expected target but no signal is observed with the antibody labeled using Zenon® technology, the signal may be increased by adding more of the Zenon® labeling reagent to increase the Fab:antibody molar ratio and/or by increasing the final concentration of the primary antibody. While a molar ratio of 3:1 is suitable in many cases, the molar ratio can be increased up to 6:1. Further increases in the molar ratio will not result in a significant increase in the signal strength.

Label loss or transfer. Because the labeled Fab fragment is not covalently coupled to the primary antibody, some loss of the labeled Fab fragments may occur over time. In applications with multiple mouse IgG primary antibodies, the Fab fragments may transfer slowly from one primary antibody to another. Samples should be subjected to an additional post-staining fixation step with formaldehyde in order to reduce Fab dissociation (see *Staining Applications*, below).

Staining Applications

The following protocols were developed at Molecular Probes for the staining of cells and tissues using primary antibodies labeled with Zenon® mouse IgG labeling reagents. The protocols may not be applicable to all experimental situations, and users are encouraged to adapt these protocols to their needs.

Immunocytochemistry

This protocol has been successfully used to stain a number of cell lines, including bovine pulmonary artery endothelial (BPAE), Muntjac and Madin-Darby canine kidney (MDCK) cells. For primary antibodies against abundant targets, forming Zenon® labeling complexes at a molar ratio of 3:1, Zenon® labeling reagent (Fab) to antibody, is suitable; for less abundant targets, a molar ratio of 6:1 is recommended.

2.1 Culture cells in suitable growth medium. Adherent cell lines intended for observation by microscopy can be grown directly on slides or coverslips; however, do not allow the slides or coverslips to dry out at any point during the staining protocol.

2.2 Fix the cell sample in 4% formaldehyde solution in PBS for 15 minutes at room temperature. Alternatively, the sample can be incubated for 15 minutes at 37°C under 5% CO₂. This latter fixation method has been observed to better preserve cytoskeletal morphology.

2.3 Wash the cell sample one or more times with PBS.

2.4 Permeabilize the cells by incubating the sample in PBS containing 0.1% Triton® X-100 for 5 minutes at room temperature. Remove the detergent solution.

2.5 Block nonspecific binding sites in the cell sample for 30 minutes in PBS containing 10% normal goat serum (NGS).

2.6 Dilute the Zenon® labeling complex (prepared according to *Zenon® Complex Formation*, above) to the desired working concentration in PBS containing 10% NGS, and apply a sufficient volume to immerse the cell sample. If staining with more than

one Zenon® labeling complex, all complexes can be combined into a single staining solution; sequential labeling is not necessary. In systems where the concentration of primary antibody has already been optimized for use with a secondary detection reagent, use the primary antibody at a concentration 1.5–3-fold higher with the Zenon® staining method. Otherwise, the appropriate working concentration will need to be determined empirically.

2.7 Incubate the cell sample for 30–60 minutes at room temperature. Zenon® labeling reagents with fluorophores should be protected from light during the incubation. When staining is complete, wash the cell sample one or more times with PBS.

2.8 Perform a second fixation of the cell sample in 4% formaldehyde solution in PBS for 15 minutes at room temperature. When the fixation is complete, wash the cell sample one or more times with PBS. Without this second fixation step, the Zenon® Fab fragment, because it is not covalently coupled to the primary antibody, can dissociate with time, and can even transfer to another primary antibody in multicolor applications. This step is optional for experiments using only a single primary antibody for labeling purposes, but signal intensity is generally superior with the additional fixation.

2.9 If desired, counterstain the cell sample with a nucleic acid stain (e.g., Hoechst 33342, DAPI or SYTOX® Green stains) or with other labeling reagents (e.g., fluorophore-labeled phalloidin), and wash with PBS after any additional staining procedure.

2.10 If necessary, deposit the cells on a slide or coverslip, or place in a flow tube. For a sample to be viewed by microscopy, mount in a suitable antifade mounting medium (e.g., use the ProLong® Gold antifade reagent; P36930). Analyze or view the sample by the desired method.

Immunohistochemistry

This protocol is designed for 10–12 µm frozen sections of tissues fixed in 4% formaldehyde via transcardial perfusion. *Do not use glutaraldehyde-containing fixatives.* The use of glutaraldehyde increases sample autofluorescence and results in excessive crosslinking of proteins, which can impede antibody access and can denature antigens. Where possible, use thinner sections to promote thorough and more rapid antibody penetration into the tissue. Zenon® labeling complexes have not been tested or optimized for sections embedded in paraffin. Note that certain low-abundance antigens may not be detected as easily using Zenon® labeling technology as they would be using conventional secondary detection methods, and the Zenon® method is not recommended for very low-abundance target molecules in tissue samples.

In parallel to the labeling of an experimental sample, a negative-control sample should be processed in an identical protocol, but with the omission of the primary antibody, in order to assess nonspecific labeling. A positive control sample labeled by conventional secondary detection methods can be included during pilot studies to confirm the specificity of the staining pattern observed using the Zenon® labeling method.

3.1 Allow the slide containing the frozen tissue section to thaw at room temperature (usually 15–30 minutes). Draw a rectangle around the section with a PAP pen, if desired.

3.2 Rehydrate the section in PBS for 15 minutes at room temperature. Do not allow the section to dry at any point during the remaining steps of this staining protocol.

3.3 Permeabilize the section in PBS containing 0.2% Triton X-100 (PBT) for 20 minutes at room temperature.

3.4 Block nonspecific binding sites with PBT containing 0.2% BSA and 5% heat-inactivated normal goat serum for 30 minutes at room temperature. Alternatively, 1% BSA can be used in place of the serum. Remove the blocking buffer once the incubation is complete.

3.5 Dilute the Zenon® labeling complex (prepared according to *Zenon® Complex Formation*, above) in PBT to the desired working concentration. In a horizontal staining dish, apply 100–250 µL of the working solution to cover the section. For initial experiments, test a range of different primary antibody concentrations and different molar ratios of Zenon® labeling reagent (Fab) to antibody. A molar ratio of 6:1 is the recommended starting point for tissue staining, and, in general, Zenon® labeling methods require a higher concentration of primary antibody than would be used in a standard secondary detection scheme. It is important to include detergent in the staining buffer in order to promote antibody penetration. Buffer PBT, as described above, contains 0.2% Triton X-100; however, the detergent choice and concentration (generally 0.1–0.5%) may need to be optimized to strike a balance between tissue permeabilization and the possibility of the antigen being sensitive to detergent extraction.

3.6 Incubate the section with the staining solution in a humidified chamber for 1–2 hours at room temperature. Zenon® labeling reagents with fluorophores should be protected from light during the incubation. If staining with more than one Zenon® labeling complex, all complexes can be combined into a single staining solution; sequential labeling is not necessary. *Do not use longer incubation times.* With longer incubation times, the Zenon® Fab fragment can dissociate from the primary antibody, resulting in poor signal-to-noise ratios and possible crosslabeling of primary antibodies in multicolor staining applications.

3.7 Remove the staining solution and wash the section in PBT for 10–15 minutes at room temperature. Repeat the wash an additional two times. The wash steps should be performed in an upright Coplin staining jar.

3.8 Wash the stained tissue section for 5 minutes in PBS at room temperature. Repeat.

3.9 Perform a second fixation of the tissue section in 4% formaldehyde solution in PBS for 15 minutes at room temperature. When the fixation is complete, wash the section one or more times with PBS. Without this second fixation step, the Zenon® Fab fragment, because it is not covalently coupled to the primary antibody, can dissociate with time, and can even transfer to another primary antibody in multicolor applications.

3.10 If desired, counterstain the tissue sample with a nucleic acid stain (e.g., Hoechst 33342, DAPI or SYTOX Green stains) or with other labeling reagents (e.g., fluorophore-labeled phalloidin), and wash with PBS after any additional staining procedure.

3.11 Mount the section in a suitable antifade mounting medium (e.g., use the ProLong® antifade reagent; P36930) and observe the fluorescence using appropriate filters.

Flow Cytometry

The following protocol was developed for the staining of peripheral blood mononuclear cells and can be readily adapted to other cell types or to whole blood samples. For primary antibodies against abundant targets, forming Zenon® labeling complexes at a molar ratio of 3:1, Zenon® labeling reagent (Fab) to antibody, is suitable; for less abundant targets the molar ratio should be increased. Note that an Fc blocking step is not required when using antibodies labeled with Zenon® technology.

4.1 Prepare the cell suspension or blood sample as preferred.

4.2 Add 100 µL of cultured cells (at 1×10^6 cells/mL) or 100 µL of a whole blood sample to a flow cytometry tube.

4.3 (Optional) Fix the cell sample in 4% formaldehyde solution in PBS for 15 minutes at room temperature. Wash the cell sample with PBS, centrifuge the sample to pellet the cells and decant the supernatant. Resuspend the cells in the residual buffer. **NOTE:** This and the following step are required only when detecting internal cell markers. If using antibodies against cell surface and internal markers in the same experiment, perform the staining of the cell surface markers first, then follow this step and steps 4.4–4.7 to stain the internal markers.

4.4 (Optional) Permeabilize the cell sample by incubating the sample in PBS containing 0.1% Triton X-100 for 5 minutes at room temperature. Wash the cell sample with PBS, centrifuge the sample to pellet the cells and decant the supernatant. Resuspend the cells in the residual buffer.

4.5 Label 1 µg of the primary antibody according to *Zenon® Complex Formation*, above.

4.6 Add the entire Zenon® labeling mixture (prepared in step 4.5) to the sample tube. **NOTE:** If using whole blood samples, the lysis of red blood cells can be performed either before or after cell staining (steps 4.6–4.7).

4.7 Incubate the sample for 30 minutes at room temperature. Zenon® labeling complexes with fluorophores should be protected from light during the incubation. When complete, wash the cells once with PBS, centrifuge the sample to pellet the cells, and decant the supernatant. Resuspend the cells in the residual buffer. A post-staining fixation step, required for imaging applications, is not required for flow cytometry applications.

4.8 Analyze the sample using appropriate instrument parameters.

Product Notes

[A] The Fab:antibody ratio is the important factor when determining the amount of Zenon® mouse IgG labeling reagent to use in the labeling protocol. In all Zenon® Mouse IgG Labeling Kits, the Zenon® labeling reagent is provided at a concentration of 200 µg/mL based upon the mass of the Fab fragment. A

Fab fragment has a molecular weight of ~50 kDa, compared to ~150 kDa for an intact IgG; thus, 5 µL of any Zenon® labeling reagent mixed with 1 µg of mouse IgG antibody produces a Fab:antibody molar ratio of 3:1.

[B] When adjusting either the amount of antibody to be labeled or the Fab:antibody molar ratio, it is important to always use equal volumes of Zenon® labeling reagent and Zenon® blocking reagent. For example, if the amount of Zenon® labeling reagent used for a reaction is increased to 10 µL, then the amount of Zenon® blocking reagent should also be increased to 10 µL. Note that adding 10 µL of the Zenon® labeling reagent for each microgram of antibody (yielding a molar ratio of 6:1) will often increase the measured signal intensities by approximately 50%. Further increases in the molar ratio tend to yield smaller increases in intensity.

[C] If the concentration of mouse IgG is not specified, contact the supplier and obtain at least an approximate IgG content. If a concentration cannot be obtained from the supplier, a titration with the Zenon® labeling reagent will need to be performed to obtain the optimum labeling conditions.

[D] Monoclonal antibodies from suppliers are generally provided as a purified IgG fraction, as ascites fluid or as hybridoma

supernatant. Primary antibodies that have not been purified can still be labeled using the Zenon® mouse IgG labeling reagents and do not require the removal of nonspecific IgGs or serum proteins. The appropriate amount of the Zenon® labeling reagent to add in step 1.2 should be determined by using the total IgG mass in the sample to be labeled; thus, 5 µL of the Zenon® mouse IgG labeling reagent should be used for each µg of IgG. Nonspecific IgGs will be labeled in addition to the specific IgG; however, the labeled nonspecific IgGs should not stain the sample appreciably and will be washed away during the staining procedure.

[E] If mouse cells or tissues are not being used, and if the application requires only a single labeled antibody, steps 1.4 and 1.5 of the labeling protocol (addition of the blocking reagent) can be omitted. The level of uncomplexed Zenon® labeling reagent is sufficiently low that it should not cause high background levels.

[F] If long-term storage of the labeling complexes is desired, stop the labeling procedure at step 1.3 of the protocol. At this stage the complexes can be stored at 4°C for several weeks with the addition of 2 mM sodium azide as a preservative. When ready to use the conjugate, complete steps 1.4 and 1.5 of the protocol before applying the conjugate to the experimental sample.

Product List Current prices may be obtained from our website or from our Customer Service Department.

Z25000	Zenon® Alexa Fluor® 350 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25100	Zenon® Alexa Fluor® 350 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25200	Zenon® Alexa Fluor® 350 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25013	Zenon® Alexa Fluor® 405 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25113	Zenon® Alexa Fluor® 405 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25213	Zenon® Alexa Fluor® 405 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25001	Zenon® Alexa Fluor® 430 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25002	Zenon® Alexa Fluor® 488 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25102	Zenon® Alexa Fluor® 488 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25202	Zenon® Alexa Fluor® 488 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25003	Zenon® Alexa Fluor® 532 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25004	Zenon® Alexa Fluor® 546 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25104	Zenon® Alexa Fluor® 546 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25204	Zenon® Alexa Fluor® 546 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25005	Zenon® Alexa Fluor® 555 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25105	Zenon® Alexa Fluor® 555 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25205	Zenon® Alexa Fluor® 555 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25006	Zenon® Alexa Fluor® 568 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25106	Zenon® Alexa Fluor® 568 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25206	Zenon® Alexa Fluor® 568 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25007	Zenon® Alexa Fluor® 594 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25107	Zenon® Alexa Fluor® 594 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25207	Zenon® Alexa Fluor® 594 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25020	Zenon® Alexa Fluor® 610-R-Phycoerythrin Mouse IgG, Labeling Kit *10 labelings*	1 kit
Z25008	Zenon® Alexa Fluor® 647 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25108	Zenon® Alexa Fluor® 647 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25208	Zenon® Alexa Fluor® 647 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25021	Zenon® Alexa Fluor® 647-R-Phycoerythrin Mouse IgG, Labeling Kit *10 labelings*	1 kit
Z25121	Zenon® Alexa Fluor® 647-R-Phycoerythrin Mouse IgG _{2b} , Labeling Kit *10 labelings*	1 kit
Z25221	Zenon® Alexa Fluor® 647-R-Phycoerythrin Mouse IgG _{2b} , Labeling Kit *10 labelings*	1 kit
Z25009	Zenon® Alexa Fluor® 660 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25109	Zenon® Alexa Fluor® 660 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25209	Zenon® Alexa Fluor® 660 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25010	Zenon® Alexa Fluor® 680 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25110	Zenon® Alexa Fluor® 680 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25210	Zenon® Alexa Fluor® 680 Mouse IgG _{2b} , Labeling Kit *50 labelings*	1 kit
Z25022	Zenon® Alexa Fluor® 680-R-Phycoerythrin Mouse IgG, Labeling Kit *10 labelings*	1 kit

Zenon® Mouse IgG Labeling Kits

Z25011	Zenon® Alexa Fluor® 700 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25030	Zenon® Alexa Fluor® 700-Allophycocyanin Mouse IgG, Labeling Kit *10 labelings*	1 kit
Z25031	Zenon® Alexa Fluor® 750-Allophycocyanin Mouse IgG, Labeling Kit *10 labelings*	1 kit
Z25050	Zenon® Alkaline Phosphatase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25150	Zenon® Alkaline Phosphatase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25250	Zenon® Alkaline Phosphatase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25051	Zenon® Allophycocyanin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25151	Zenon® Allophycocyanin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25251	Zenon® Allophycocyanin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25052	Zenon® Biotin-XX Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25152	Zenon® Biotin-XX Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25252	Zenon® Biotin-XX Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25053	Zenon® DSB-X™ Biotin Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25042	Zenon® Fluorescein Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25054	Zenon® Horseradish Peroxidase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25154	Zenon® Horseradish Peroxidase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25254	Zenon® Horseradish Peroxidase Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25040	Zenon® Marina Blue® Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25043	Zenon® Oregon Green® 488 Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25041	Zenon® Pacific Blue™ Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25156	Zenon® Pacific Blue™ Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25055	Zenon® R-Phycoerythrin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25155	Zenon® R-Phycoerythrin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25255	Zenon® R-Phycoerythrin Mouse IgG, Labeling Kit *25 labelings*	1 kit
Z25045	Zenon® Texas Red®-X Mouse IgG, Labeling Kit *50 labelings*	1 kit
Z25060	Zenon® Tricolor Mouse IgG, Labeling Kit #1 *for green, orange and deep red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25070	Zenon® Tricolor Mouse IgG, Labeling Kit #2 *for blue, green and red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25080	Zenon® Tricolor Mouse IgG, Labeling Kit #3 *for flow cytometry, 488 nm excitation* *3 x 10 labelings*	1 kit
Z25160	Zenon® Tricolor Mouse IgG, Labeling Kit #1 *for green, orange and deep red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25170	Zenon® Tricolor Mouse IgG, Labeling Kit #2 *for blue, green and red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25180	Zenon® Tricolor Mouse IgG, Labeling Kit #3 *for flow cytometry, 488 nm excitation* *3 x 10 labelings*	1 kit
Z25260	Zenon® Tricolor Mouse IgG, Labeling Kit #1 *for green, orange and deep red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25270	Zenon® Tricolor Mouse IgG, Labeling Kit #2 *for blue, green and red fluorescence imaging* *3 x 10 labelings*	1 kit
Z25280	Zenon® Tricolor Mouse IgG, Labeling Kit #3 *for flow cytometry, 488 nm excitation* *3 x 10 labelings*	1 kit

Zenon® Mouse IgG Labeling Kits

Contact Information

Further information on Molecular Probes products, including product bibliographies, is available from your local distributor or directly from Molecular Probes. Customers in Europe, Africa and the Middle East should contact our office in Paisley, United Kingdom. All others should contact our Technical Service Department in Eugene, Oregon.

Please visit our website—probes.invitrogen.com—for the most up-to-date information.

Molecular Probes, Inc.

29851 Willow Creek Road, Eugene, OR 97402
Phone: (541) 465-8300 • Fax: (541) 335-0504

Customer Service: 6:00 am to 4:30 pm (Pacific Time)

Phone: (541) 335-0338 • Fax: (541) 335-0305 • probesorder@invitrogen.com

Toll-Free Ordering for USA:

Order Phone: (800) 438-2209 • Order Fax: (800) 438-0228

Technical Service: 8:00 am to 4:00 pm (Pacific Time)

Phone: (541) 335-0353 • Toll-Free (800) 438-2209
Fax: (541) 335-0238 • probestech@invitrogen.com

Invitrogen European Headquarters

Invitrogen, Ltd.
3 Fountain Drive

Inchinnan Business Park

Paisley PA4 9RF, UK

Phone: +44 (0) 141 814 6100 • Fax: +44 (0) 141 814 6260

Email: eurolinfo@invitrogen.com

Technical Services: eurotech@invitrogen.com

Molecular Probes products are high-quality reagents and materials intended for research purposes only. These products must be used by, or directly under the supervision of, a technically qualified individual experienced in handling potentially hazardous chemicals. Please read the Material Safety Data Sheet provided for each product; other regulatory considerations may apply.

Limited Use Label License

For research use only. Not intended for any animal or human therapeutic or diagnostic use. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any Commercial Purpose, and that such collaborator agrees in writing (a) to not transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Invitrogen Corporation will not assert a claim against the buyer of infringement of the above patents based upon the manufacture, use or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture of such product. If the purchaser is not willing to accept the limitations of this limited use statement, Invitrogen is willing to accept return of the product with a full refund. For information on purchasing a license to this product for purposes other than research, contact Molecular Probes, Inc., Business Development, 29851 Willow Creek Road, Eugene, OR 97402. Tel: (541) 465-8300. Fax: (541) 335-0504.

Several Molecular Probes products and product applications are covered by U.S. and foreign patents and patents pending. All names containing the designation * are registered with the U.S. Patent and Trademark Office.

Copyright 2005, Molecular Probes, Inc. All rights reserved. This information is subject to change without notice.